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TRATAMENT FACIAL | LIFTING | REJUVENARE

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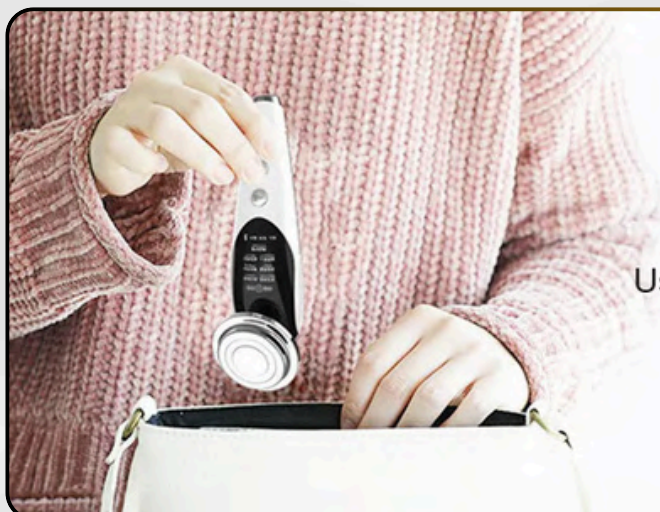
7-IN-1 PROFESSIONAL BEAUTY DEVICE

imatrend®

FACIAL TREATMENT | LIFTING | REJUVENATION

EMS & RF Technology | Photonic LED Therapy |
Needle-Free Mesotherapy | Sonic Vibration Massage

Lifting & Firming | Anti-Aging Facial Treatment |
Targeted Eye Care | Deep Pore Cleansing

USER MANUAL

TABLE OF CONTENTS

| | |
|--|-----------|
| 1. Safety Precautions | |
|  | 21 |
| 2. Product Specifications | |
|  | 22 |
| 3. Product Components | |
|  | 23 |
| 5. List of accessories | |
| 6. Technology principles | |
|  | 24 |
| 6. Operating Modes | |
|  | 25 |
| 7. Instructions for the 7 operating modes | |
|  | 26 |
| 8. Step-by-step guide for the 7 operating modes | |
|  | 27 |
| 9. Skincare Tips (Routine) | |
|  | 32 |
| 10. Additional warnings and precautions | |
|  | 33 |
| 11. Product maintenance and cleaning | |
|  | 34 |
| 12. Troubleshooting | |
| 13. Contact | |
|  | 35 |
| 14. Warranty Certificate | |
|  | 36 |

1. SAFETY PRECAUTIONS

Safety Operation Manual

- Please read this manual carefully before using the product. Keep this manual in a safe place. If you transfer the product to a third party, ensure this manual is included with the device.
- Any operation that does not comply with the instructions in this manual may lead to serious personal injury, property damage, or damage to the product itself.
- Our company shall not assume any legal liability for losses caused by failure to operate the product in accordance with the provided instructions.
- Please pay close attention to the personal safety instructions to avoid accidents or potential damage to the device.
- Exercise caution while using the product to avoid potential hazards.
- Store the product safely and keep it out of reach of children or persons unable to assess environmental risks.
- Do not store the product in high places to avoid the risk of accidental falling.

Caution:

- Do not use the device if you have had recent facial plastic surgery, or if you have subcutaneous metal implants or a pacemaker.
- Prohibited for individuals with skin diseases, allergies to metals or cosmetics, and patients with tuberculosis, hypertension, or acute medical conditions.
- Not recommended during pregnancy, breastfeeding, menstruation, or for children under 12 years of age.
- Do not use on open wounds, mucous membranes, the Adam's apple, or directly on the eyeball.
- Remove all metal jewelry, glasses, and contact lenses prior to use.
- The recommended usage time should not exceed 15 minutes per session.

2. PRODUCT SPECIFICATIONS

| Parameter | Technical Details |
|-------------------------|---|
| Product Name | 7-in-1 Professional Anti-aging and Rejuvenation Device |
| Input Voltage | 5V / 1A (USB Charging) |
| Rated Power | 4W |
| Battery Capacity | 750 mAh (Lithium-Ion) |
| Charging Time | Approx. 3 hours |
| Standby Autonomy | Approx. 45 days |
| Intensity Levels | 3 Levels (Low / Medium / High) |
| Integrated Technologies | EMS, RF, Sonic Vibration, LED Therapy (Red, Blue, Orange) |
| Treatment Areas | Face, Eyes, Lips, Neck |
| Usage | Dry skin, Wet skin |
| Skin Type | Oily, Combination, Blemish-prone, Sensitive, Normal |
| Dimensions | 160 x 50 x 54 mm (6.29 x 1.96 x 2.12 in) |
| Net Weight | 100 g (3.20 oz) |
| Package Contents | Main device, Fixing ring, Base, Charging cable, User manual |

3. PRODUCT COMPONENTS

Metal Treatment Head

Control Interface and Display

To operate the device correctly, please refer to the LCD indicators and the control button functions described below:

A. Battery Symbol Indicator

The energy status is indicated through color codes and illumination patterns:

- White: Normal battery level.
- Red: Low battery (device requires charging).
- Flashing: The device is connected and currently charging.
- Solid Light (Fixed): Charging is complete.

Full charge required before first use (approx. 3 hours).

B. BUTTON FUNCTIONS

The device features two primary buttons that control all integrated functions:

1. Intensity Button:

- **Short Press:** Adjusts the electrical current intensity levels (3 stages: **L - Low, M - Medium, H - High**).
- **Long Press:** Activates or Deactivates the Vibration function. This feature is versatile and can be toggled for each operating mode according to your personal comfort.

Intensity Button

2. ON/OFF & MODE Button:

- **Long Press (3 seconds):** Powers the device ON or OFF.
- **Short Press:** Cycles through the 7 operating modes displayed on the screen.

ON/OFF & MODE Button

4. LIST OF ACCESSORIES

- Main Device
- Fixing Ring (pre-mounted on the device)
- Charging / Support Base
- USB Data Cable
- User Manual

5. TECHNOLOGY PRINCIPLES

The Ima Trend Anti-aging and Facial Rejuvenation Device is an advanced beauty tool that integrates EMS microcurrents (galvanic ions), thermal massage vibrations, and photonic therapy (Red, Blue, and Orange light).

This technological synergy is engineered to achieve deep pore cleansing, visible tissue rejuvenation, and a long-lasting facial lifting effect.



Before beginning your treatment, it is important to understand the advanced technologies working for your skin's health:

1. **EMS Microcurrent (Galvanic Ions):** Utilizes the principle of gentle impulse currents (anode and cathode) to effectively attract, soften, and discharge deep-seated oil and impurities from the pores.
2. **Vibrations:** Employs physical vibrations and precise electrical stimulation to increase skin elasticity and boost cellular metabolism.
3. **Thermal Massage:** Promotes optimal absorption of skincare products, helps smooth skin texture, and maintains a firm, elastic complexion.
4. **LED Photonic Therapy:**
 - **Red Light (640nm-650nm):** Penetrates deep into the dermis, stimulating blood circulation and collagen production to reduce wrinkles and firm the skin. It encourages cellular regeneration.
 - **Blue Light (415nm):** Provides rapid anti-inflammatory effects. It reduces sebum secretion and eliminates acne-causing bacteria, significantly improving oily or problematic skin conditions.
 - **Orange Light (600nm):** Reaches 1-2 mm into the skin to effectively treat red spots, soothe inflammation, and promote metabolism for a brighter, more radiant complexion.






6. OPERATING MODES

This device offers customized solutions for various skin concerns by combining heat, vibration, and LED therapy:

- **CLEAN:** The warming sensation helps open pores for deep cleansing. The result is refined pores, diminished post-acne spots, and a revitalized look for tired or dull skin.
- **RF LED:** Utilizes red light to activate cellular elasticity. It promotes protein and collagen synthesis, working in synergy with thermal massage and vibrations to leave the skin smooth and firm.

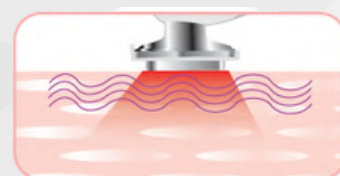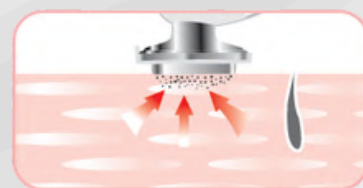- **ACNE-T:** Blue light is ideal for targeting acne, as it effectively inhibits skin inflammation and excessive oil (sebum) secretion, visibly improving pore size and skin texture.
- **WHITEN (BRIGHTENING):** Orange light restores cellular energy and accelerates skin metabolism. It revitalizes the skin to restore its natural glow and achieve a brighter, more radiant complexion.
- **IMPORT (MESOTHERAPY):** Microcurrents and thermal action drive skincare products (serums, essences) deep into the skin. Combined with LED therapy and vibrations, this mode ensures that the lower layers of the dermis absorb nutrients, leaving the skin deeply hydrated and protected.
- **EMS UP (LIFTING):** EMS microcurrent massage boosts blood circulation while enhancing skin elasticity and firmness. This mode is designed to lift facial contours and leave your skin in its peak condition.
- **EYE CARE (ÎNGRIJIREA OCHILOR):** A thermal vibration massage specifically engineered for the delicate eye area. It relieves fatigue, drives eye cream deep into the tissue, and promotes maximum absorption for a refreshed and luminous look.

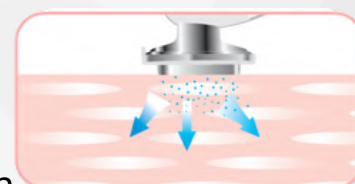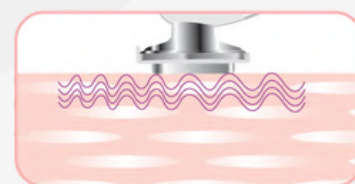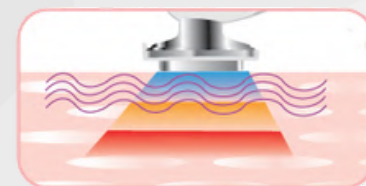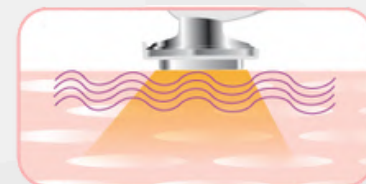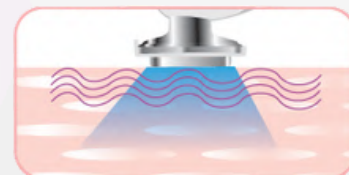7. INSTRUCTIONS FOR THE 7 OPERATING MODES

General Safety and Efficiency Rules:

1. **Preparation:** Fully charge the device before its first use.
2. **Sensitivity Test:** Test the device on a small area of your skin (such as the jawline) to check for sensitivity. Always start with Intensity Level 1.
3. **The Human Circuit:** For EMS and Ion technology to function, you must keep your fingers in contact with the Ion Sensor Strips (side rails) while the metal treatment head touches your skin. Always ensure your skin remains moist during treatment.
4. **Correct Cleansing:** In CLEAN mode, use only toner or water-based makeup remover. Do not use oil-based cleansers or greasy makeup removal lotions for exfoliation.

5. **Sensor Contact:** If using the device on another person, you must maintain physical contact (e.g., touch their skin with your other hand) to complete the electrical circuit and activate the ion/EMS functions.
6. **Auto-Shutoff:** For your safety and to conserve battery, the device will automatically power off after 5 minutes of inactivity in the same mode.


IMPORTANT NOTE:

BEFORE USE, PLEASE CAREFULLY READ SECTION NO. 10: IMPORTANT WARNINGS AND SAFETY MEASURES ON PAGE 33.


8. STEP-BY-STEP GUIDE FOR THE 7 OPERATING MODES

[CLEAN] Mode (Deep Cleansing)

This mode must be used with the fixing ring and a cotton pad.

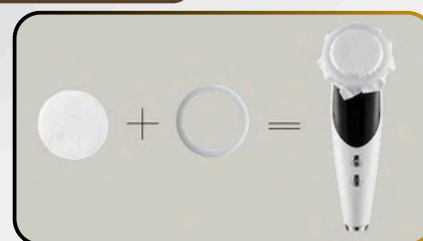- **Preparation:** Place a clean cotton pad over the metal treatment head and secure it with the plastic fixing ring. Press gently until it clicks.
- **Application:** Pour makeup remover or cleansing toner onto the cotton pad until it is well-soaked (do not use oils).
- **Power On:** Long-press the Power button (3 sec.) and select the CLEAN mode.
- **Operation:** Ensure your fingers are touching the side sensor strips during use. Move the device over your face in circular motions, avoiding the eye area.

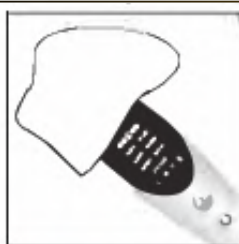1. Place a clean cotton pad

2. Fix the ring


3. Pour makeup remover or cleansing toner

4. Power On: Long-press the Power button (3 sec.) and select the CLEAN mode.


5. Ensure your fingers are touching the side sensor strips during use.

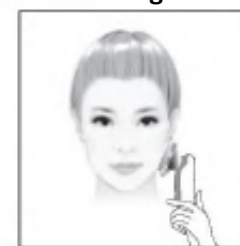6. Clean the face

Pro Tips for CLEAN Mode:

- **DO NOT** use oil-based cleansers: These can clog pores or damage the fixing ring. Use only water-based products (micellar water, toner).
- **Direction of movement:** Always move the device from the bottom up and from the center of the face outward.
- **Frequency:** We recommend using this mode 2-3 times per week for optimal results without compromising the skin's natural barrier.

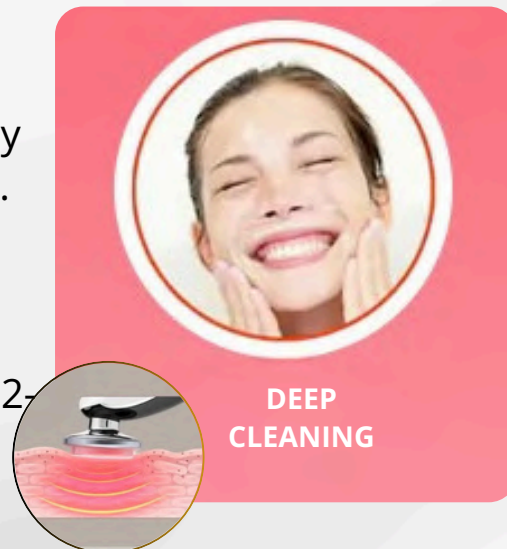[RF LED / ACNE-T / WHITEN] Modes

These modes are used with toner, serum, or lotion **WITHOUT** the fixing ring and **WITHOUT** a cotton pad. They rely on direct contact between the metal treatment head and hydrated skin.

- **Application:** Evenly spread your favorite toner, serum, or lotion over your face (ensure the skin is thoroughly moist).
- **Power On:** Long-press the Power button and select the desired mode based on your skin's needs (Red/Blue/Orange).
- **Operation:** Grip the device so your hand is in constant contact with the side sensor strips (ion sensors).
- **Ascent Massage:** Massage the face using slow, upward motions to promote the lifting and regeneration effects.

Specific Details for Each Mode:

- **RF LED (Red Light):** Ideal for collagen stimulation. Best used in combination with anti-aging serums.

- **ACNE-T (Blue Light):** Ideal for calming irritated skin. Best used with water-based products or lightweight anti-acne treatments.
- **WHITEN (Orange Light):** Ideal for evening out the skin tone. Best used with Vitamin C serums or brightening products.

IMPORTANT NOTE: Do not look directly into the LED light during use. If you feel the skin has become dry during the massage, apply more product to maintain optimal conductivity.

3. [IMPORT] Mode (Mesotherapy / Absorption)

This mode is used WITHOUT the fixing ring and WITHOUT a cotton pad. It is specifically engineered to work perfectly with sheet masks.

- **Application:** Apply a sheet mask to your face (alternatively, use a generous layer of gel-mask or essence).
- **Power On:** Long-press the ON/OFF button (3 sec.) and select the IMPORT mode using short presses.
- **Operation:** Massage your face over the sheet mask or directly on the skin using upward motions, focusing on areas that require deep hydration. The microcurrent will drive nutrients deep into the dermis.
- **Final Step:** Remove the sheet mask after the 5-minute session and massage any remaining essence into the skin.


1. **Apply facial mask on the face.**

2. **Power On:** Long-press the Power button (3 sec.) and select the IMPORT mode.

3. **Hold the device, touching the side sensor strips and massage the face from bottom to top.**

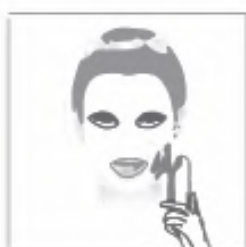4. **Remove the facial mask and clean the face.**

Pro Tips for IMPORT Mode:


- The "Lock-in" Effect: This mode uses positive and negative ions to "push" hydrating molecules through the skin barrier. It is 10 times more effective than applying cream with your fingers.
- Ideal for Sheet Masks: If using a sheet mask, move the device over it until you feel the mask's fabric beginning to dry. This indicates that the serum has been successfully transferred into the skin.
- Temperature Matters: The thermal massage in this mode helps relax the pores, facilitating the "import" of nutrients. If the cosmetic product dries too quickly, mist your face with a little toner or thermal water.
- Correct Order: Use this mode AFTER cleansing (CLEAN Mode) and BEFORE your final heavy night cream to "seal" the vitamins inside.
- No Waste: You can also use this mode for the neck and décolletage—areas that often reveal age and require a collagen and hydration "boost."

4. [EMS UP] Mode (Lifting)

This mode uses microcurrents to stimulate facial muscles. It is used WITHOUT the ring or cotton pad. Ideal for firmness and redefining facial contours.

- **Application:** Evenly apply a lifting cream, day/night cream, conductive gel, or hydrating lotion to your face.
- **Power On:** Long-press the ON/OFF button (3 sec.) and select the EMS UP mode using short presses.
- **Sensor Contact:** Ensure your fingers are touching the side sensor strips during use to activate the microcurrents.
- **Operation:** Always massage the face from the bottom up, along the jawline and over the cheekbones using upward motions. The microcurrents will stimulate the facial muscles for a visible lifting effect.

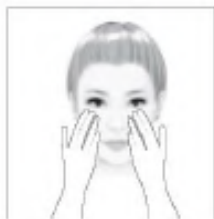1. Apply toner/essence/
lotion on the face evenly

2. Power On: Long-press the
Power button (3 sec.) and
select EMS mode.

3. Hold the device, touching
the side sensor strips
and massage the face
from bottom to top.

Pro Tips for EMS Mode:

- The "Tingling" Sensation: It is completely normal to feel slight pricking or tingling. This indicates that the microcurrents are effectively working on the muscles. If the sensation is too intense, reduce the intensity to the L (Low) level.
- The 45° Rule: Hold the metal treatment head at an angle so it gently "hooks" the jawline, lifting it towards the ear.
- Focus on Critical Areas: Use slow motions along the facial "V-line" (jaw), cheekbones, and neck (avoiding the central neck area/Adam's apple).
- Hydration is Key: If you feel no sensation, your skin may be too dry. Apply more serum or lotion, as microcurrents require a moist environment to conduct effectively.
- Consistency: For a visible and long-lasting lifting effect, use this mode consistently (3-4 times per week).

5. [EYE CARE] Mode

Specially designed for the sensitive under-eye area and eyelids. Use with a specialized eye cream or periocular serum.

- **Application:** Gently tap eye cream into the skin around the eyes.
- **Power On:** Long-press the ON/OFF button and select EYE CARE mode with short presses.
- **Operation:** Keep your fingers on the sensors and gently massage from the inner corner to the outer corner of the eyes, applying light pressure.
- **Note:** The microcurrent for this mode is disabled by default for safety. If you require more intense stimulation, you must manually activate it by short-pressing the Intensity button while in this mode.

1. Apply eye cream or eye essence around the eye evenly


2. Power On: Long-press the Power button (3 sec.) and select EYE CARE mode.


3. Hold the device, and make sure the hand is in touch with the sensor, massage from inner corners to outer corners of the eye area and gently pat and press into skin.

EYE CARE EMS micro-current is off by default, and can be turned on by pressing the Intensity button.

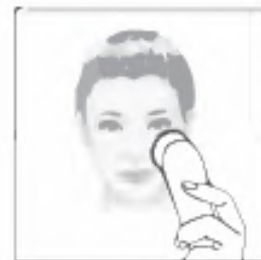Pro Tips for EYE CARE Mode:

- **Manual Microcurrent Activation:**
By default, the microcurrent is disabled to protect the sensitive eye area. If you desire a more intense eyelid lifting effect, short-press the Intensity button to activate it manually.
- **No Excessive Pressure:** Do not press the device directly against the eyeball! Movements should be light and focused on lymphatic drainage along the orbital bone (the area under the eye and above the eyebrow).
- **No Dark Circles:** The device's warmth helps dilate blood vessels and promotes lymphatic drainage—ideal for reducing morning puffiness.
- **Optimal Temperature:** The thermal massage helps your eye cream penetrate deep into the skin layers. You will notice your skin becoming significantly brighter and more hydrated.
- **Direction Matters:** Always massage toward the temples. This movement helps "lift" the gaze and eliminate fluids accumulated under the eyes.


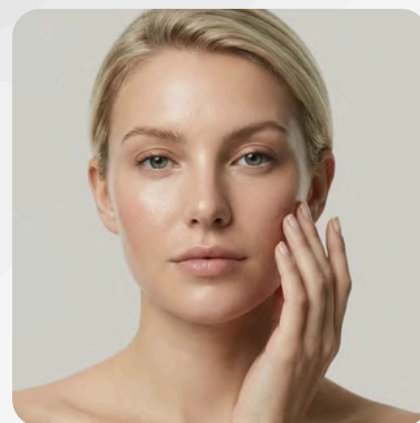
9. SKINCARE TIPS (THE IMA TREND ROUTINE)

Morning Treatment (Makeup Prep):

- **Cleansing:** Wash your face with water or a gentle cleansing gel.
- **Preparation:** Apply your toner or daytime serum.
- **Treatment:** Select RF LED / ACNE-T / WHITEN mode (choose the light based on your skin's needs: Red for toning, Blue for calming, Orange for radiance).

- **Finishing:** Apply your day cream and makeup (your skin will be significantly more receptive and smooth).

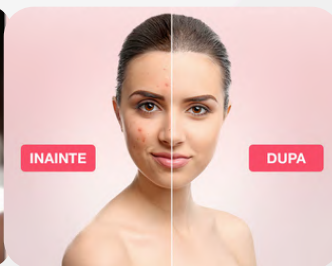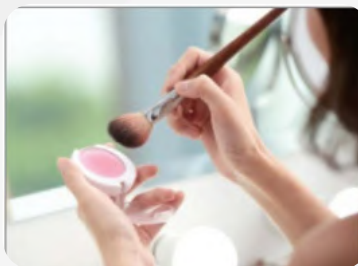Evening Treatment (Vitality Restoration):

- **Step 1:** Deep Cleansing ([CLEAN] Mode): Use the fixing ring and a cotton pad with makeup remover or micellar water. (Recommended 2–3 times per week to avoid compromising the skin's natural barrier).
- **Step 2:** Photonic Treatment ([RF LED / ACNE-T / WHITEN] Mode): Apply toner or essence. The light therapy assists in cellular regeneration throughout the night.
- **Step 3:** Intense Hydration ([IMPORT] Mode): Apply a sheet mask. Use the device over the mask to drive nutrients deep into the dermis.
- **Step 4:** Recontouring ([EMS UP] Mode): After removing the mask, use the remaining serum on your face to perform the lifting massage.
- **Step 5:** Specific Care ([EYE CARE] Mode): Apply eye cream and massage gently to drain dark circles before sleep.

10. ADDITIONAL WARNINGS AND PRECAUTIONS

SAFETY AND CONTRAINDICATIONS

- DO NOT use the device if you fall into any of the following categories:
- Procedures and Implants: Individuals with recent facial plastic surgery, subcutaneous metal implants, pacemakers, or silicone/plastic parts embedded in the treatment area.
- Medical Conditions: Patients with hypertension, infectious diseases, tuberculosis, acute conditions, mental illnesses, or atopic skin diseases.
- Health Status: Individuals in poor health, under the influence of alcohol, or taking sleeping aids.
- Sensitivities: Severe allergies to metals or specific cosmetics (please consult a physician before use).
- Special Categories: Women who are pregnant, breastfeeding, or during their menstrual cycle. This product is NOT intended for individuals under the age of 12.

PROHIBITED AREAS for Use:

- Open wounds, cuts, surgical scars, inflamed or purulent (infected) areas.
- Mucous membranes (oral cavity, genital area), the eyeballs, or directly over the Adam's apple.
- Areas with viral warts or recently waxed/shaved skin.

IMPORTANT PRE-OPERATING WARNINGS

- Preparation: You must remove all metal jewelry (earrings, necklaces, watches), glasses, and contact lenses. Metal contact with the device may cause accidental burns or microcurrent malfunction.
- Eye Protection: The LED light can be intense. Keep your eyes closed when working near the face and DO NOT look directly into the LED light source.
- Treatment Head: The metal surface is precision-finished for delicate contact and can be easily scratched if handled carelessly. Fine scratches occurring during use do not affect performance and are not covered by the warranty.
- Usage Time: Do not exceed 15 minutes per treatment session.
- Proper Use: Do not use the metal head to pry open or force other parts of the device.
- Battery Maintenance: If the device is not used for an extended period, charge it fully at least once every 6 months to protect the lifespan of the lithium battery.

11. PRODUCT MAINTENANCE AND CLEANING

- To ensure optimal operation, always turn off the device before cleaning.
- DO NOT WASH: The device is not waterproof. Do not submerge it or rinse it under running water.
- Cleaning the Body: Wipe the device after each use with a soft, dry cloth or a very slightly damp cloth.
- Cleaning the Fixing Ring: Remove the ring and wash it separately with water when dirty. Dry it completely before reattaching it to the device.
- WARNING Regarding Corrosive Substances: Do not use alcohol, benzene, nail polish remover, or other corrosive solutions. These can cause parts to crack, discolor the casing, or lead to permanent sensor failure.

- **Battery** : If the device remains unused, you must charge it at least once every 6 months to keep the lithium cells active and maintain their health.
- **Recycling**: At the end of its lifespan, this device must be disposed of in accordance with local regulations for electronic waste (WEEE) and lithium batteries. Do not dispose of it with regular household waste.

NOTE: If you have any concerns regarding your health status, please consult a physician before the first use.

12. TROUBLESHOOTING

| Issue | Possible Cause | Solution |
|---|--|---|
| Device will not power on | The battery is completely discharged. | Charge the device for at least 3 hours. |
| No vibration or microcurrent sensation | Not touching the side sensors or the skin is too dry. | Ensure your hand is in contact with the side sensor strips and apply more hydrating product to your face. |
| Burning or stinging sensation | Intensity level is too high or insufficient product on the skin. | Reduce the intensity level to (L) and apply more gel or serum. |
| Device turns off by itself | The 5-minute Auto-Shutoff safety function. | This is normal. Restart the device if you wish to continue the treatment. |

CONTACT

If you encounter any issues or have suggestions, please feel free to contact us directly via our customer support phone line.

+4 0755. 622 .628
(WhatsApp, Viber)

Instagram:
@imatrend.ro

office@imatrend.ro

Mon - Fri:
09:00 AM - 06:00 PM

https://imatrend.ro

WARRANTY CERTIFICATE

COMMERCIAL WARRANTY FOR BUSINESS CUSTOMERS: 12 MONTHS
CONFORMITY WARRANTY FOR INDIVIDUALS: 24 MONTHS

The manufacturing, assembly, and performance of this product have been tested and meet the required quality standards.
The product has been approved for release.

Name

Address

.....

Product: **IMA Trend 7-in-1 Professional Beauty Device**

Purchase Date

Phone

The warranty does not apply in the following situations: unauthorized repairs or modifications, incorrect use, negligence, abuse, collision, or liquid ingress, use of non-original or incompatible accessories, breaking, modifying, or removing labels (including authenticity labels), exceeding the warranty period, lack of a valid document confirming the purchase date, damage caused by human factors or unavoidable forces (natural disasters), and normal wear and tear of the product.

TRAVEL-FRIENDLY

Easy to carry—enjoy your skincare routine anytime, anywhere.

imatrend[®]